



Green University of Bangladesh

Department of Computer Science and Engineering

Lab report-4

Course Title: Chemistry Laboratory

Course code: CHE-102

Date of Submission: 27.12.2020

Submitted to:	Submitted by :
Name : MR.FORKAN SAROAR	Name : Nazifa Alam Nowrin
Designation : Lecturer	ID : 193002103
Department : Textile	Department : CSE
Green university of Bangladesh	Green university of Bangladesh

Experiment number: 04

Experiment name: Standardization of Sodium Thiosulphate solution with standard Potassium Dichromate solution.

Objectives:

1. To study the strength of thiosulphate.
2. To study oxidation reduction titration.
3. To study the liberation and its titration.

Learning Outcome:

After completing this experiment the students will be able to:

1. Determine the strength of thiosulphate with the standard dichromate solution & liberated iodine.
2. Observe the end point by color change with starch indicator

Theory:




Standardization is the process by which the strength of a solution is determined with the help of a standard solution. A solution of known concentration is called a

standard solution. This experiment is done by means of titration. In presence of a suitable indicator, a chemical substance that detects the end point of reaction by changing its color, the volumetric analysis in which a standard solution is added in

another solution (whose strength is unknown) to reach its end point to determine the strength of that solution is called titration.

Titration involving iodine or dealing with liberated iodine in chemical reaction is called iodimetry and iodometry respectively. This reaction is iodometric because

iodine is obtained from KI. The reactions of this experiment are:

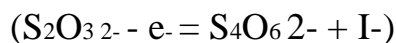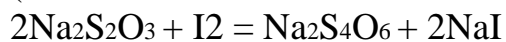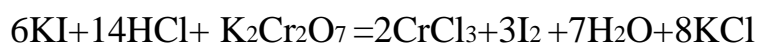


Here $\text{K}_2\text{Cr}_2\text{O}_7$ is an oxidizing agent and I^- is a reducing agent. Again in the second

reaction I_2 is an oxidizing agent and $\text{S}_2\text{O}_3^{2-}$ is a reducing agent.

In the 2nd Step of the reaction a specific indicator is used that is “Starch”- which has a significant effect on iodine.

$\text{Starch} + \text{I}_2 = \text{Starch-Iodine}$ (blue colour)

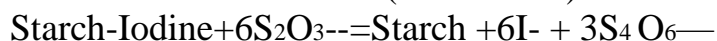

If to a solution containing a little iodine, some starch solution is added and $\text{Na}_2\text{S}_2\text{O}_3$ is run in from the burette, the blue color of the starch-iodine complex will disappear from the solution as soon as all the iodine has been reduced to iodide ion.

Apparatus:

1. Conical flask,
2. Burette,
3. Pipette,
4. Volumetric flask,
5. Stand,
6. Funnel

Chemicals:

1. $\text{Na}_2\text{S}_2\text{O}_3$,
2. $\text{K}_2\text{Cr}_2\text{O}_7$,
3. KI,
4. NaHCO_3
5. HCl (Concentrated),
6. Starch (Indicator)

Data and calculation:

(Standardization of $\text{Na}_2\text{S}_2\text{O}_3$ solution with standard $\text{K}_2\text{Cr}_2\text{O}_7$ solution)

Number of Objects	Volume of $\text{K}_2\text{Cr}_2\text{O}_7$ (ml)	Burette (ml)	reading	Volume of $\text{Na}_2\text{S}_2\text{O}_3$ (ml)	Average volume of $\text{Na}_2\text{S}_2\text{O}_3$ (ml)
Initial Reading	Final Reading				
1	10	0	5	5	
2	10	5	10.2	5.2	5.1
3	10	10.2	15.3	15.3	

Calculation:

We know,

$$6V_{\text{red}} \times S_{\text{red}} = V_{\text{ox}} \times S_{\text{ox}}$$

Here,

$$V_{\text{K}_2\text{Cr}_2\text{O}_7} = 10\text{ml}, S_{\text{K}_2\text{Cr}_2\text{O}_7} = 0.5\text{M}$$

$$V_{\text{Na}_2\text{S}_2\text{O}_4} = 5.1\text{ml}, S_{\text{Na}_2\text{S}_2\text{O}_4} = ?$$

$$S_{\text{Na}_2\text{S}_2\text{O}_4} = 6(10 \times 0.5) / 5.1$$

Result:

Determined strength of $\text{Na}_2\text{S}_2\text{O}_3$ solution is = 5.88M

Discussion:

As the color change of the titration of $\text{Na}_2\text{S}_2\text{O}_3$ with $\text{K}_2\text{Cr}_2\text{O}_7$ is very confusing,
the end point of the titration may not have been properly determined. This may
be
the cause of error